

Fresh Air Hood for Autopsy Area ZU 10/900

The fresh air hood has especially been developed for the area above the Autopsy table. Stable construction made of stainless steel for installation in the ceiling or independent installation. Optimal and draught free work is guaranteed by adjustable air-conduction. A regular airflow-emission prevents the ascent of pollutants. The fresh air hood is to be connected to the side-mounted ventilation system. The integrated special-lamps (4x58 Watt) are providing regular light, similar to daylight, with colour-light temperature of approx. 5500°K.

Technical data and dimensions, see back page



Fresh Air Hood build into the ceiling




Autopsy Table with Fresh Air Hood

Subject to technical modifications as required by product improvements.

Fresh Air Hood for Autopsy Area ZU 10/900

Dimensions W/D/H:	1085 x 2400 x 400 mm
Air capacity:	approx. 750 m ³ /h
Loss of pressure:	approx. 40 pa
Air-conduction:	adjustable
Lighting:	4 x 58 Watt
Brightness:	approx. 2000 Lux at a distance of 1,5m
Temperatur of colour:	approx. 5500° Kelvin
Voltage:	230 Volt / 50 Hz
Exhaust duct diameter:	2 x 200 mm
Material:	Stainless steel Material no. 1.4301 - ANSI 304
Type:	ZU 10/900
Article no:	010.010.690

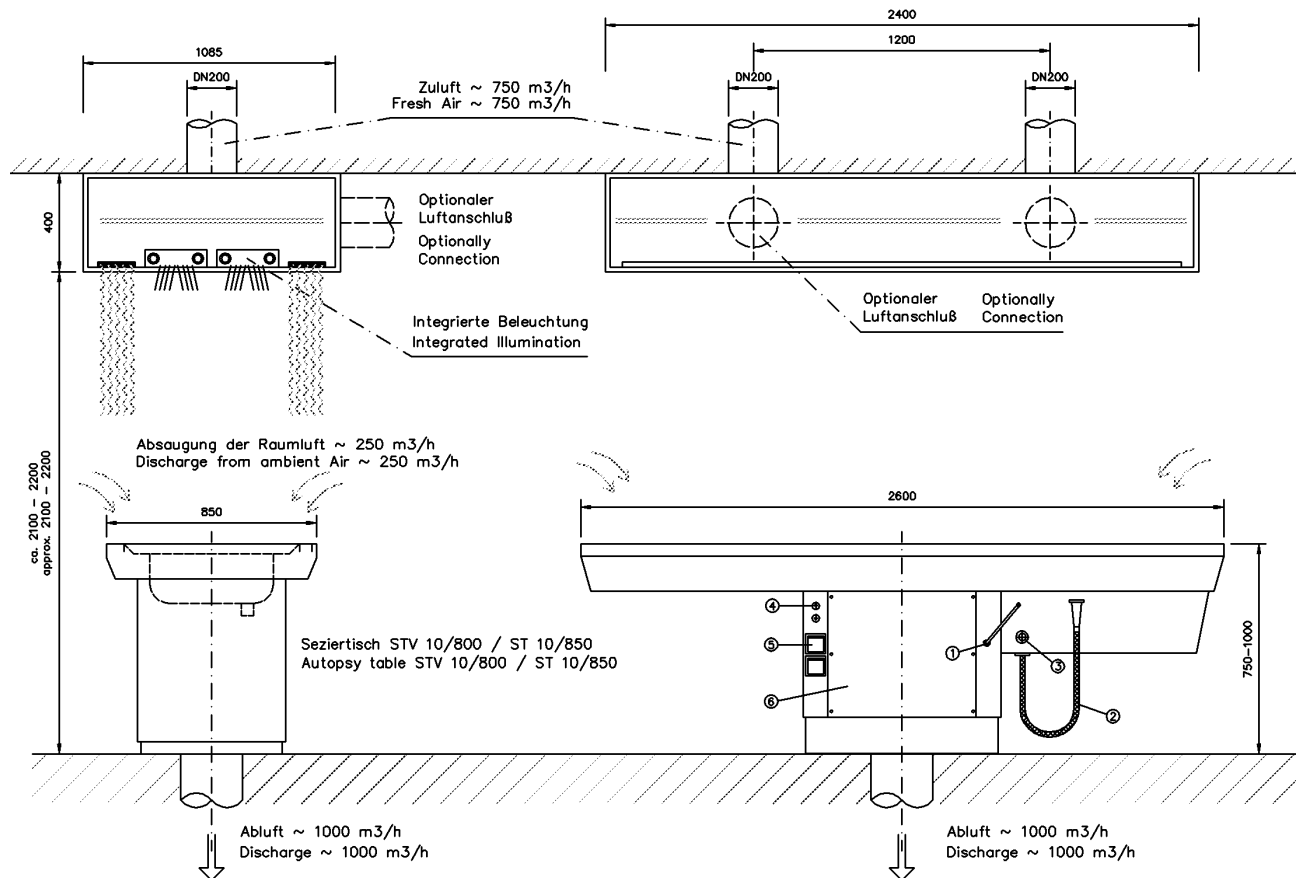

Subject to technical modifications as required by product improvements.